

# Embgenix™ PGT-A Kit: Preparation of Whole Genome Amplification

---

## Author

[Opentrons](#)

## Categories

---

- NGS Library Prep
  - Embgenix™ PGT-A

## Description

---

This protocol performs the [Embgenix™ PGT-A Kit: Preparation of Whole Genome Amplification](#), parts V-VI:B.


## Modules

- [Opentrons Temperature Module \(GEN2\)](#)







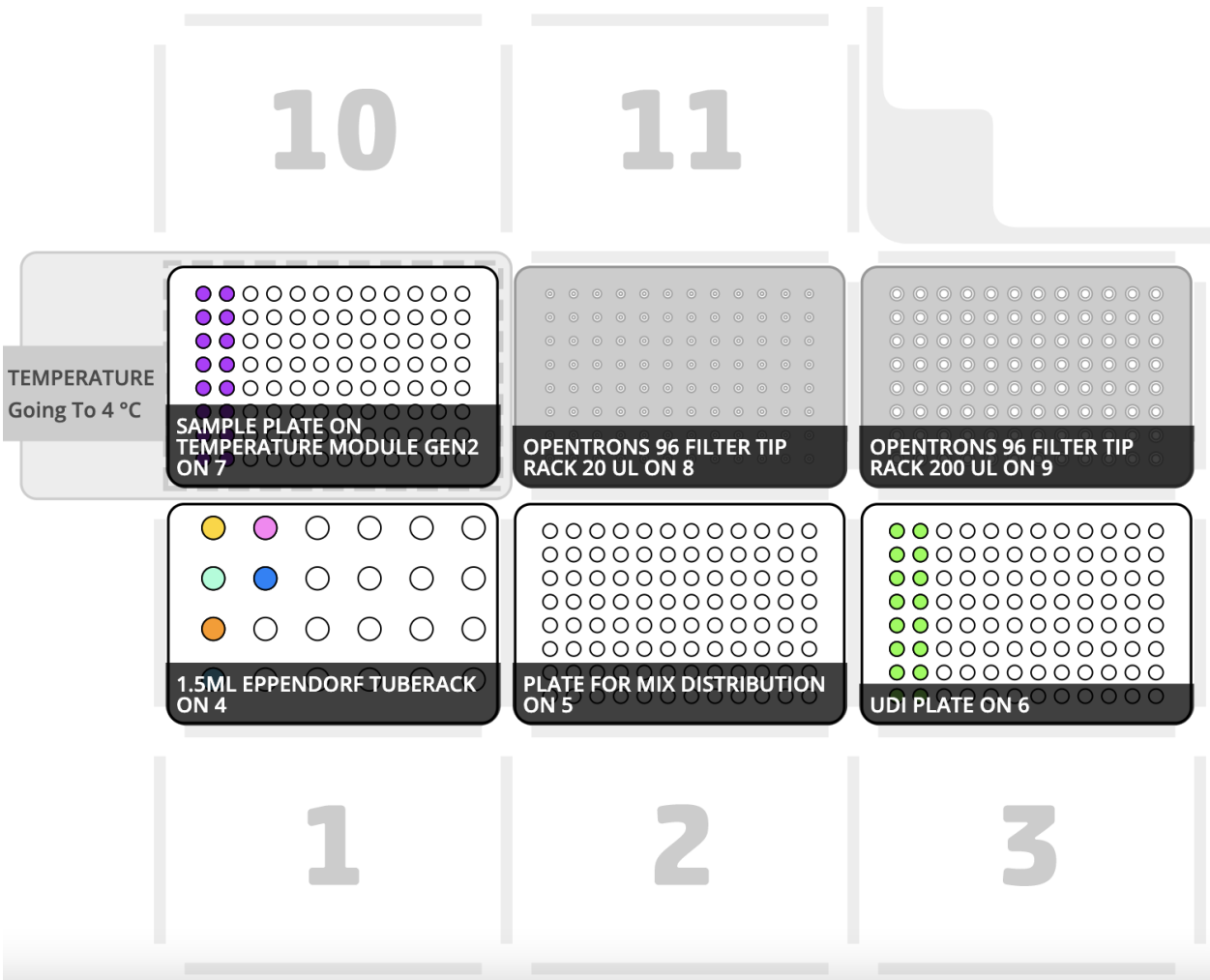


## Labware

- [Opentrons 96 Well Aluminum Block with NEST Well Plate 100 µL](#)
- [Opentrons 24 Tube Rack with Eppendorf 1.5 mL Safe-Lock Snapcap](#)

Deck Setup



## Reagent Setup

 SAMPLE
 CELL EXTRACTION MASTERMIX
 WHOLE GENOME AMPLIFICATION MASTERMIX
 WD1
 WD2
 LIBRARY PREP MASTERMIX
 LIBRARY AMPLIFICATION MASTERMIX
 STEM LOOP ADAPTERS
 UDI

- Columns 1-2 on Temperature Module Plate (slot 7):  
Starting sample
    - A1 on Tuberack (slot 4): Cell Extraction Master Mix
    - B1 on Tuberack (slot 4): Whole Genome Amplification Master Mix
  - C1 on Tuberack (slot 4): WD2
  - D1 on Tuberack (slot 4): WD2
  - A2 on Tuberack (slot 4): Library Prep Master Mix
  - B2 on Tuberack (slot 4): Library Amplification Master Mix
- Mix
- Columns 5-6 on Mix Distribution Plate (slot 5): Stem Loop Adapters
- Loop Adapters
- Columns 1-2 on UDI Plate (slot 6): UDI

## Process

1. Input your protocol parameters above.
2. Download your protocol and unzip if needed.
3. Upload your custom labware to the [OT App](#) by navigating to `More > Custom Labware > Add Labware`, and selecting your labware files (.json extensions) if needed.
4. Upload your protocol file (.py extension) to the [OT App](#) in the `Protocol` tab.
5. Set up your deck according to the deck map.
6. Calibrate your labware, tiprack and pipette using the OT App. For calibration tips, check out our [support articles](#).
7. Hit "Run".

## Additional Notes

If you have any questions about this protocol, please contact the Protocol Development Team by filling out the [Troubleshooting Survey](#).