

Cell Culture

Author

[Opentrons](#)

Categories

- Sterile Workflows
 - Cell Culture

Description

This protocol performs a custom 12-factor cell culture assay in a 96-deepwell plate.


Labware

- [Opentrons 6 Tube Rack with Falcon 50 mL Conical](#)
- [USA Scientific 96 Deep Well Plate 2.4 mL #1896-2000](#)
- Opentrons 96 Filter Tip Rack 200 μ L
- Opentrons 96 Filter Tip Rack 1000 μ L







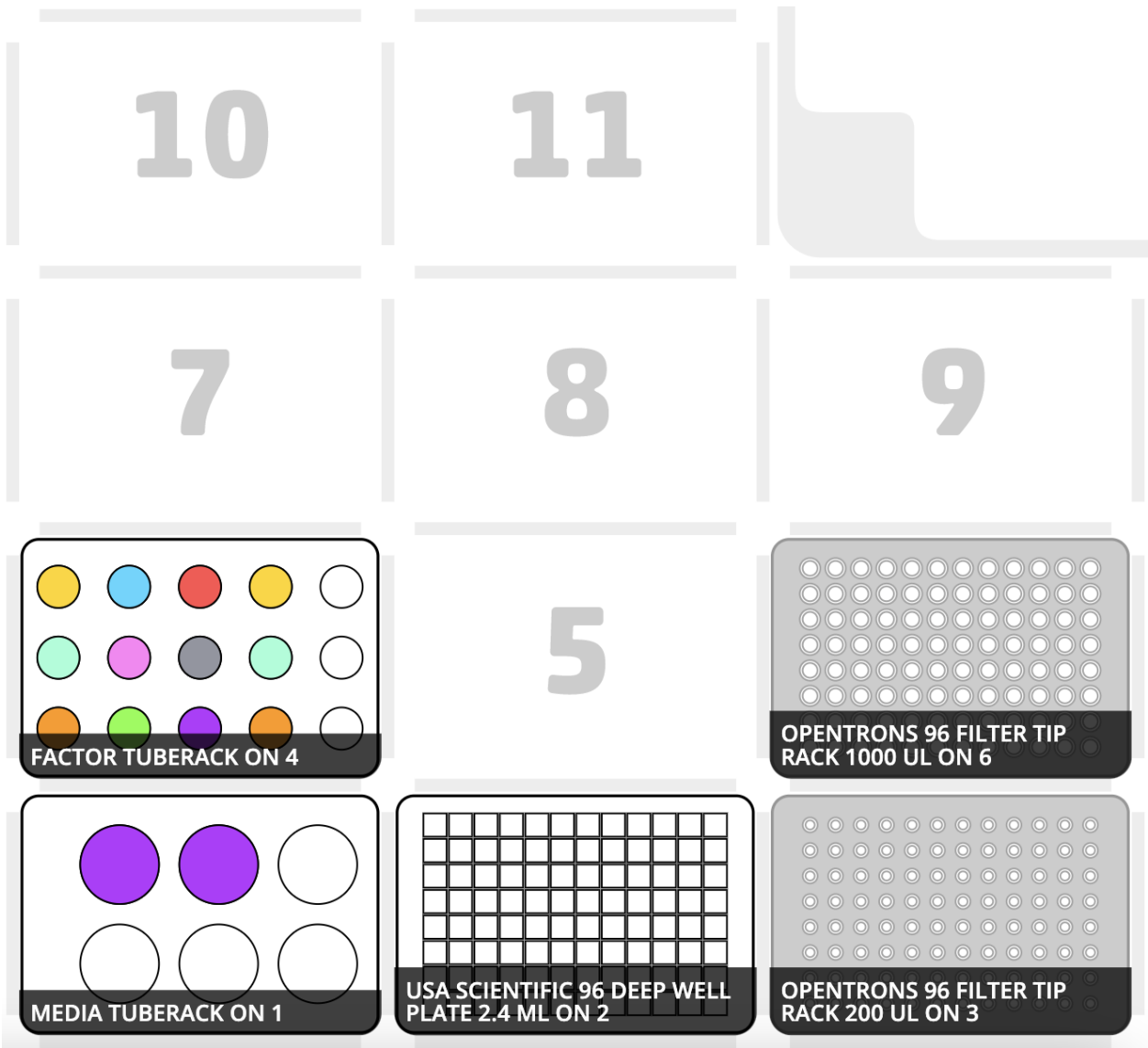






Pipettes

- [Opentrons P1000 Single Channel Electronic Pipette \(GEN2\)](#)
- [Opentrons P300 Single Channel Electronic Pipette \(GEN2\)](#)

Deck Setup



Reagent Setup

Liquids	
	MEDIA
	FGF8
	BGJ
	BMP4/8
	LDN
	WNT3A
	IWP2
	R-SPONDIN
	GDNF
	SPP86
	LARA
	AGN193109
	GAMMA XX

Process

1. Input your protocol parameters above.
2. Download your protocol and unzip if needed.
3. Upload your custom labware to the [OT App](#) by navigating to `More > Custom Labware > Add Labware`, and selecting your labware files (.json extensions) if needed.
4. Upload your protocol file (.py extension) to the [OT App](#) in the `Protocol` tab.
5. Set up your deck according to the deck map.
6. Calibrate your labware, tiprack and pipette using the OT App. For calibration tips, check out our [support articles](#).
7. Hit "Run".

Additional Notes

If you have any questions about this protocol, please contact the Protocol Development Team by filling out the [Troubleshooting Survey](#).